



**NÜVE SANAYİ MALZEMELERİ İMALAT VE TİCARET A.Ş.**

**NF 800 / 800 R**

**MULTI-PURPOSE BENCH  
TOP CENTRIFUGES**

**USER'S MANUAL**



Z14. K25 254 Rev. No: 07 Rev.Date: 01/2016

**MANUFACTURER :**

**NÜVE SANAYİ MALZEMELERİ  
İMALAT VE TİCARET A.Ş.**

Saracalar Mah. Saracalar Kümeevleri No: 4/2

Akyurt 06750 Ankara-TURKEY

TEL: +(90) 312 399 28 30 (pbx)

FAX: +(90) 312 399 21 97

E-mail: sales@nuve.com.tr

**WARRANTY CERTIFICATE**

1. Nüve warrants that the equipment delivered is free from defects during material and workmanship. This warranty is provided for a period of two years. The warranty period begins from the delivery date.
  2. Warranty does not apply to parts normally consumed during operation or general maintenance or any adjustments described in the operating instructions provided with the equipment.
  3. Nüve does not accept any liability in the case where the goods are not used in accordance with their proper intent.
  4. The warranty may not be claimed for damages occurred during the shipment, for damages resulting from improper handling or use, the defects in maintenance, negligence, bad functioning of auxiliary equipment, in the case of force majeure or accident and incorrect power supply.
  5. In the event of failure, Nüve shall be under no liability for any injury, or any loss or damage as the result of the failure other than the guarantee conditions.
- **BEFORE OPERATING THE INSTRUMENT THIS MANUAL SHOULD BE READ CAREFULLY.**
  - **THE VALIDITY OF THE GUARANTEE IS SUBJECT TO THE OBSERVATION OF THE INSTRUCTIONS AND PRECAUTIONS DESCRIBED IN THIS MANUAL.**
  - **INFORMATION CONTAINED IN THIS DOCUMENT IS THE PROPERTY OF NÜVE. IT MAY NOT BE DUPLICATED OR DISTRIBUTED WITHOUT PERMISSION.**

Dear Nüve User,

We would like to take this opportunity to thank you for preferring this Nüve product. Please read the operating instructions carefully and keep them handy for future reference.

Please detain the packing material until you see that the unit is in good condition and it is operating properly. If an external or internal damage is observed, contact the transportation company immediately and report the damage. According to ICC regulations, this responsibility belongs to the customer.

While you are operating the instrument please;

- obey all the warning labels,
- do not remove the warning labels,
- do not operate damaged instrument,
- do not operate the instrument with a damaged cable,
- do not move the instrument during operation.

In case of a problem contact your Nüve agent for an authorized service or maintenance.

The validity of the guarantee is subject to compliance with the instructions and precautions described in this manual.

Nüve reserves the right to improve or change the design of its products without any obligation to modify previously manufactured products.

Information contained in this document is the property of Nüve. It may not be duplicated or distributed without its permission.

**PLEASE REGISTER ONLINE TO VALIDATE WARRANTY:**

To register your warranty online, please visit our web page [www.nuve.com.tr](http://www.nuve.com.tr) and fill in **WARRANTY REGISTRATION FORM**.

## TABLE OF CONTENTS

<b>SECTION 1 - INTRODUCTION .....</b>	<b>4</b>
1.1 USE AND FUNCTION.....	4
<b>SECTION 2 - TECHNICAL SPECIFICATIONS .....</b>	<b>5</b>
2.1 TECHNICAL SPECIFICATIONS TABLE .....	5
2.2 ACCESSORIES FOR NF 400 / 400R .....	5
2.3 GENERAL PRESENTATION .....	6
2.4 ROTOR SELECTION TABLE.....	7
2.5 PRECAUTIONS AND USAGE LIMITATIONS .....	7
<b>SECTION 3 - SYMBOLS.....</b>	<b>8</b>
<b>SECTION 4 - INSTALLATION .....</b>	<b>8</b>
4.1 LIFTING AND TRANSPORT .....	8
4.2 UNPACKING .....	9
4.3 POSITIONING .....	9
4.4 MAINS SUPPLY .....	9
<b>SECTION 5 - OPERATING UNIT .....</b>	<b>10</b>
5.1 OPERATING .....	10
5.2 DISPLAYS AND CONTROL PANEL .....	10
5.2.1 Introduction to Control Display Panel Components and Functions.....	10
5.3 MAKING A NEW PROGRAM .....	12
5.4 RUNNING THE PROGRAM.....	13
5.4.1 SELECTING THE MEMORIZED PROGRAM.....	13
5.5 SAFETY INTERLOCK SYSTEM .....	14
5.6 IMBALANCE DETECTION SYSTEM .....	14
5.7 MANUAL LID OPENING .....	14
<b>SECTION 6 - OPERATING PRINCIPLES .....</b>	<b>15</b>
6.1 PREPARATION OF THE ROTOR TO RUN.....	15
6.2 LOADING .....	15
6.3 ROTOR INSTALLATION.....	16
6.4 DRIVE SYSTEM .....	17
<b>SECTION 7 - CLEANING AND PERIODIC MAINTENANCE .....</b>	<b>17</b>
7.1 PERIODIC MAINTENANCE.....	17
7.2 STERILIZATION .....	17
7.3 CORROSION INFORMATION .....	17
7.3.1 Chemical Corrosion.....	18
7.3.2 Stress Corrosion.....	18
7.4 CLEANING .....	18
7.5 ELECTRICITY .....	19
<b>SECTION 8 - DISPOSAL MANAGEMENT CONCEPT .....</b>	<b>19</b>
<b>SECTION 9 - TROUBLESHOOTING .....</b>	<b>19</b>
<b>SECTION 10 - ELECTRICAL CIRCUIT DIAGRAMS.....</b>	<b>21</b>
10.1 ELEKTRICAL CIRCUIT DIAGRAM OF NF 800 .....	21
10.2 ELECTRICAL CIRCUIT DIAGRAM OF NF 800 R.....	22
<b>SECTION 11 - WARNING LABEL .....</b>	<b>23</b>

# SECTION 1

## INTRODUCTION

### 1.1 USE AND FUNCTION

The NF 800 / NF 800 R centrifuges are designed to separate liquids of different densities by applying centrifugal force.

Different types of rotors, swing-out and angle rotors, can be fitted to the centrifuges.

The microprocessor control system of the centrifuges enables the user to observe and to control the speed (RPM or RCF), the acceleration and braking rates, run time, program number and temperature of the bowl for NF 800 R. The microprocessor control system also displays the RPM or RCF value during the run.

The ventilation system of the NF 800 limits temperature in the bowl between the ambient temperature + 3 °C and ambient temperature + 12 °C depending on the rotor type. The NF 800 R is equipped with a cooling system which provides to operate between -9 °C and + 40 ° C.

**The NF 800/800R centrifuges are manufactured according to the following standards, EN 61010-1, EN 61010-2-020, EN 61000-6-3, EN 50419.**

**This device is in compliance with WEEE Regulation.**

## SECTION 2

### TECHNICAL SPECIFICATIONS

#### 2.1 TECHNICAL SPECIFICATIONS TABLE

	NF 800	NF 800 R
Max Speed	RA 200 rotor: 4.100 rpm, RS24 rotor: 14.000 rpm	
Max. RCF	RA 200 rotor: 3.045 x g, RS24 rotor : 18.405 x g	
Max. Capacity	RA 200 rotor: 4x 200 ml, RS 24 rotor : 24 x 1.5/2 ml	
Control System	Microprocessor controlled	
Display	High Visibility Digital	
No. of Memories	10	
Speed Range	500-14000 rpm	
Speed Step	10 rpm	
Speed Accuracy	± 20 rpm	
Timer Range	01-99 min + Hold Position	
Timer Step	1 min	
Acceleration Rate	0: Slowest 9: Fastest	
Braking Rate	0: Slowest 9: Fastest	
Power supply	230 V, 50 Hz	
Temp. range		-9° C / 40 ° C
Refrigerant Liquid		R134a
Temp. Accuracy		1°C
Temp. Performance	Limited to ambient temperature +12° C at max speed	+4°C at max. speed
Max. Power	500 W	800 W
Ext. Dimension	455 x 545 x 400	770 x 555 x 400
Packing Dim.	590 x 520 x 490	870 x 650 x 550
Net / Gross Weight	42 / 46 Kg	75 / 80 Kg

#### 2.2 ACCESSORIES FOR NF 800 / 800R

- **A 08 062** Puller for rotor removal
- **B 50 013** RA 200 swing- out rotor 4\*200 ml , 4.100rpm , 3.045 x g
- **G 51 026** Set of 4 buckes for RA 200
- **G 51 033** Set of 4 buckets with sealing cap
- **G 51 032** Set of 4 200 ml PP bottles with cap
- **G 51 031** Set of 4 insert 1x100 ml
- **G 51 030** Set of 4 inserts 1x50ml conical
- **G 51 028** Set of 4 inserts 7x15 ml
- **G 51 029** Set of 4 inserts 3x15 ml conical
- **G 51 027** Set of inserts 7x5/7 ml
- **B 50 015** RS 450 Angle rotor 30x15ml 4.100 rpm 3.007xg
- **G 03 012** Adaptor for 1.5 / 2 ml tubes  
(Should be ordered according to the capacity of the rotor)
- **G 03 014** Adaptor for 13 x 100 mm vacuumed / non - vacuumed rotor  
(Should be ordered according to the capacity of the rotor)

- **G 03 015** Adaptor for 13 x 75 mm vacuumed / non- vacuumed rotor (Should be ordered according to the capacity of the rotor)
- **B 50 014** RS 24 Angle rotor 24 x1,5/2 ml 14.000 rpm , 15.405xg
- **A 14 005** Adaptor 500/800 µl
- **A 14 006** Adaptor 200 µl PCR
- **A 14 007** Adaptor 250/400/700 µl
- **B 50 017** RS 300 Angle rotor 6x50 ml , 9.000 rpm , 9.418xg
- **G 51 024** Set of 6 adaptors 1x15 ml cocinal
- **G 51 025** Set of 6 rubber base adaptors for 50 ml cocinal tubes
- **B 50 021** MP 200 Microtitle 4.100 rpm , 2.349xg

### 2.3 GENERAL PRESENTATION



Figure 1

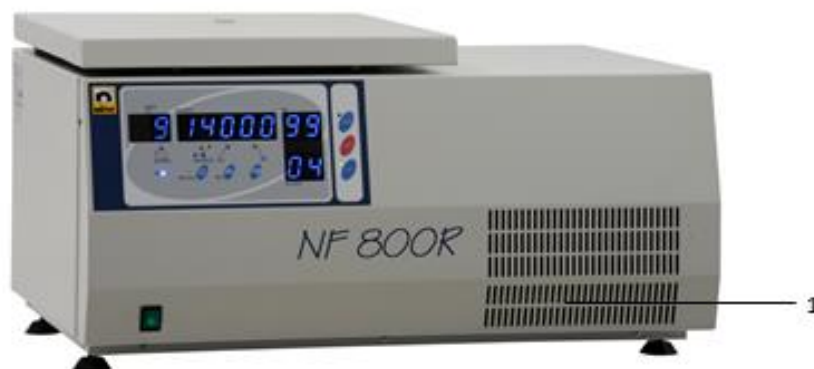


Figure 2

1	LID	7	On / off switch
2	Lock pins	8	Bowl
3	Gas spring	9	Electrical cables inlet and fuses
4	Gasket	10	Manual LID opening hole
5	Rotor	11	Cooling unit
6	Display and control panel		

## 2.4 ROTOR SELECTION TABLE

Rotor	Description	Capacity	Max. Tube Dim. Mm.	Rad mm	Max. Speed Rpm	Max RCF xg
RA 200	Swing-Out Rotor	4x 200 ml	Ø 57	162	4.100	3.044
	Set of 4 buckets	200 ml	Ø 57	162	4.100	3.044
	4 PP bottles with cap	200 ml	Ø 56,5	162	4.100	3.044
	Set of 4 inserts	1x100 ml	Ø 45	162	4.100	3.044
	Set of 4 inserts	1x50 ml Conical	Ø 30	162	4.100	3.044
	Set of 4 inserts	7x15 ml	Ø 17	162	4.100	3.044
	Set of 4 inserts	3x15 ml Conical	Ø 17	162	4.100	3.044
	Set of 4 inserts	7x5/7	Ø 13,5	162	4.100	3.007
MP 200	Microtitre Rotor	2 x 3 microtitre		125	4100	2349
RS 300	Angle Rotor	6 x 50 ml.	29.5	104	9000	9418
	Set of 6 adapters	1 x 15 ml. conical	17			
	Set of 6 rubber base adapters	1 x 50 ml. conical	29.5			
RS 450	Angle rotor	30x15 ml	Ø 17	160	4.100	3.007
	Adaptor	1 x 5 ml	Ø 13	160	4.100	3.007
	Adaptor	1 x 7 ml	Ø 13	160	4.100	3.007
	Adaptor	1 x 1.5/2 ml	Ø 11	160	4.100	3.007
RS 24	Angle Rotor	24x 1,5 / 2 ml	Ø 11	84	14.000	18.405
	Set of 24 adapters	1x500 / 800µl	Ø 8	84	14.000	18.405
	Set of 24 adapters	1x 200 µl PRC	Ø 6.5	84	14.000	18.405
	Set of 24 adapters	1x250/400/700µl	Ø 6.5	84	14.000	18.405

**Note:** Insert set is not needed for 200 ml tubes.

## 2.5 PRECAUTIONS AND USAGE LIMITATIONS

- Do not use the device for any purpose other than the usage purpose.
- Prior to first use, the user's manual should be read and the device is only to be used by authorized and trained personnel. Only authorized technical personnel handle the product in case of any failure.
- The working bench should be durable to the device weight and vibration isolated.
- Ensure that the rotor is placed correctly prior to usage.
- According to the standard IEC 61010-2-020, anyone and any hazardous materials should not be in the 300 mm safety zone while centrifuge is running.
- Do not move the device while it is running.
- Do not open the lid while rotor is spinning.
- Apply the manual lid opening procedure in the case of power cut or in the case of any error.





- Use only the the spare parts, rotors and accessories which are supplied by NUVE.
- Load the rotor according to the explanations in the user's manual.
- Start the device after ensuring the rotor is loaded correctly.
- Do not use the centrifuge in areas which are in explosive danger.
- Do not centrifuge the explosive, flammable, radioactive, corrosive materials and the materials which may react with each other.
- The centrifuge and the rotor are not microbiologically leak-proof. Use tubes with leak-proof covers, if hazardous, toxic and pathogenic microorganisms are centrifuged.
- Do not use corrosive materials which may be harmful for the device integrity, rotor and accessories.
- Do not use rotors and accessories with corrosion and mechanical damages.
- Mains supply should be appropriate to power of the device and grounded.
- Use tubes whose sizes are suitable to the rotor and accessories.
- Tubes which are used in the centrifuge should not be deformed by the effect of the centrifuge force.
- Use glass tubes to balance, if glass tubes are used. Use plastics tubes to balance, if plastics tubes are used.
- Do not start the device unless tubes are in balance.
- Imbalance loading may cause mixing the samples, broken tubes, and damages on the rotor and motor shaft.

**IF MENTIONED WARNINGS ARE NOT CONSIDERED, NÜVE WILL NOT BE RESPONSIBLE FROM THEIR RESULTS.**

## SECTION 3

### SYMBOLS

	<p>Symbol in the operating instructions: Attention, general hazard area. This symbol refers to safety relevant warnings and indicates possibly dangerous situations. The non-adherence to these warnings can lead to material damage and injury to personal.</p>
	<p>Symbol in the operating instructions: This symbol refers to important circumstances.</p>

## SECTION 4

### INSTALLATION

#### 4.1 LIFTING AND TRANSPORT

Because of the heavy weight of the centrifuge, all lifting and transport must be carried out using proper handling equipment. The centrifuge should be lifted from underneath and never be turned over.

## 4.2 UNPACKING

Remove the packing cardboard box and the second nylon packing around the centrifuge. The below written are provided with the instrument, please check them;

- 1 ea. user's manual
- 1ea. Warranty certification
- 1 ea. power cable
- 1 ea. manual lid opening tool
- 1 ea. 24 socket wrench
- 1 ea. Wrench
- Lubricant oil

## 4.3 POSITIONING

- Check that no damage has occurred during transport.
- Check that the positioning is suitable for users.
- Lift the centrifuge (use proper handling equipment if necessary) underneath and carry it to its place. Check that the centrifuge is stable on its four pads.
- The bench-top must be rigid enough to support the weight and vibrations.
- Leave 30 cm. free space on the every side of the centrifuge.
- Open the lid by using unlocking tool and check no substances left in the bowl. (See the manual lid opening section 5.7 )
- Make sure that the centrifuge does not occupy the utilization space of nearby equipment or do damage to them.
- The user should be able to follow the centrifuge during the operation.



According to the standard IEC 61010-2-020, anyone and any hazardous materials should not be in the 300 mm safety zone while centrifuge is running.

### ATTENTION !!!

The centrifuges are designed to operate safely under the following conditions:

- Indoor use only
- Ambient temperature: 5 °C to 25 °C.
- Maximum relative humidity of 80% for temperature up to 22 °C.
- Maximum altitude: 2000 m.
- Temperature for maximum performance: 15 °C / 25 °C.

## 4.4 MAINS SUPPLY

- The centrifuge requires 230 V, 50 / 60 Hz.
- Please make sure that the supplied mains match the required power ratings. If no, provide an extra line to support.



**Always plug the sterilizer/ovens to properly earthed sockets.**



**A supply fitted with a circuit breaker should be used for protection against indirect contact in case of an insulation fault.**

# SECTION 5

## OPERATING UNIT

### 5.1 OPERATING

- Open the lid and check that there is no sample in the instrument.
- Ensure that the centrifuge is plug-in to correctly grounded sockets.
- Turn the centrifuges on by using On/Off switch.
- Learn the functions of the control panel (See Part 5.2.1.).
- Set the values and start the operation (See Part 5.3)

### 5.2 DISPLAYS AND CONTROL PANEL

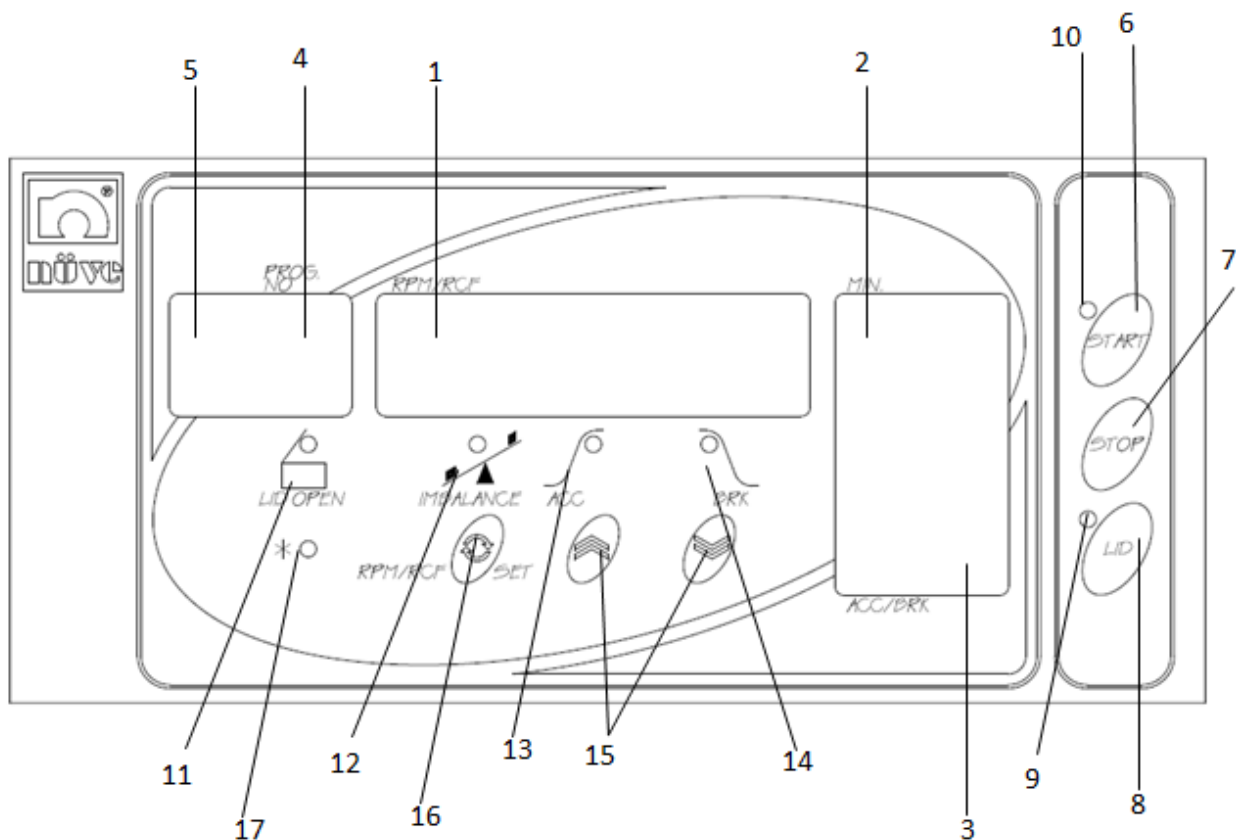


Figure 3

#### 5.2.1 Introduction to Control Display Panel Components and Functions

##### 01- Speed Display

- Display of speed (500-14.000 rpm) or display of RCF (0-18.405 x g) during the run and rotor type during programming.

##### 02- Time Display

- Display of the run time from 1 min. to 99 min. and hold position (H).

##### 03- Acceleration/Brake/Temperature Display

- Display of the acceleration and brake rates from 0 to 10, for NF 800 R it also displays the bowl

temperature during programming and the run.

**04- Program No Display**

- Display of the program number from 0 to 9.

**05- Running indicator leds**

- The 4 leds light up consecutively during the run to indicate the spinning of the rotor.

**06- START key**

- Starts the centrifugation run.

**07-STOP key**

- Stops the run by starting the braking phase manually.

**08- LID key**

- Opens the lid if the lid indicator lights up.

**09- LID Indicator Led**

- This led is on if the opening of the lid is possible.

**10- START Indicator Led**

- This led turns on as the start button is pushed and turns off as the “Lid Indicator”(9) turns on after the rotor has stopped spinning.

**11-LID OPEN Indicator Led**

- It turns on if the lid remains open or is not properly closed.

**12- Imbalance Indicator Led**

- It is activated if imbalance occurs. Brake is applied immediately.

**13- Acceleration Indicator Led**

- It turns on while the acceleration rate is being programmed and while the rotor is accelerating.

**14- Braking Indicator Led**

- This indicator turns on while the braking rate is being programmed and while the rotor is braking.

**15-Value Increase/Decrease Keys**

- They increase or decrease values during programming stage and to select the program number during stand-by position. By pushing value increase key it is possible to see the set values during centrifugation.

**16- Function key**

- It is pushed to see the RPM or RCF value on the speed display during the run, to pass among the displays and to store the programmed values.

**17- Cooling Indicator Led\***

- It indicates that the cooling system functions. During programming, it indicates that the temperature value is being set on Acc/Br/Temp Display (3).


\* It is not valid for NF 800 model.

### 5.3 MAKING A NEW PROGRAM


Determine the rotor type, the speed and time, the acceleration and braking rates.

- Push the function key (16) once and see that the value on the program no display (4) flashes. Enter the program number by pushing the value increase and value decrease keys (15).
- Push the function key (16) again and see that the value on the speed display (1) flashes. Select the proper rotor type by pushing value increase and value decrease keys (15).

ROTOR TYPE	DISPLAY		MAXIMUM SPEED (RPM)	MAXIMUM SPEED (RCF)
RA 200 Swing Out	S	out	4.100	3.044
MP 200			4.100	2.349
RS 450 Angle	30	15	4.100	3.006
RS 300 Angle	6	50	9.000	9.418
RS 24 Angle	24	1.5	14.000	18.405

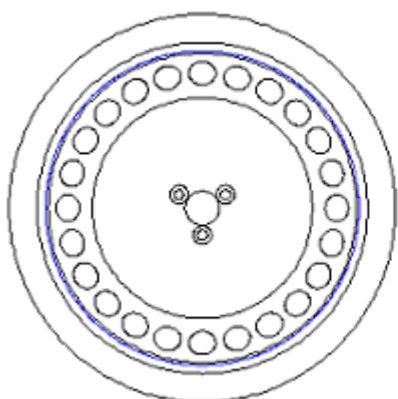

Angle Rotor

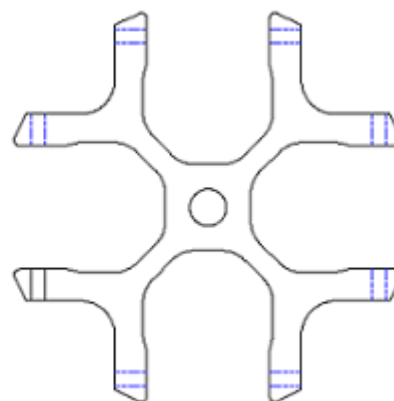

Swing-out Rotor

**NOTE :** The correct rotor type should be selected in order to see the correct RCF values during the run.

- Push the function key (16) and see that the RPM or RCF option on the speed display (1) flashes. Choose the RPM or RCF option by pushing value increase and value decrease keys (15).
- Push the function key (16) and see that the value on the speed display (1) flashes. Enter the run speed or RCF by pushing value increase and value decrease keys (15).
- Push the function key (16) and see that the value on the time display (2) flashes. Enter the run time by pushing the value increase and decrease keys (16). (H means hold position)
- Push the function key (16) and see that the acceleration led turns on and the value on the temperature display (3) flashes. Choose the acceleration value among 0,1,2,3,4,5,6,7,8 or 9 by pushing the value increase and decrease keys (15). (0: slowest 10: fastest rates)

- Push the function key (16) and see that the break led turns on and the value on the acc/break/temperature display (3) flashes. Choose the break value among 0,1,2,3,4,5,6,7,8 or 9 by pushing value increase and decrease keys.(0: slowest 9: fastest rates)
- Push the function button (16) and see that the cooling indicator led (17) turns on and the value on the acc/break/temperature display (3) flashes. Enter the bowl temperature value. (This step is not valid for NF 800 model.)
- Push the function key (16), the entered values will be saved.

You may make 10 (0-9) different programs.

## 5.4 RUNNING THE PROGRAM

- Switch on the centrifuge, see that the on/off switch led turns on.
- See that the control panel activates.
- Load your samples to the rotor and close the lid.

You can run the centrifuge in two different ways,

- You can make a new program (see section 4.3)
- You can select a program from the memory

### PRE-COOLING

While operating the NF 800 R centrifuge at sub-ambient temperatures, in order to get the maximum performance, without loading the rotor run the centrifuge at 1.500 rpm, at a temperature which is 3° C lower than the required temperature and at hold position until you see that the temperature inside the bowl drops to the required temperature. Then, stop the centrifuge, load your samples and start centrifugation.

#### 5.4.1 SELECTING THE MEMORIZED PROGRAM

Select the requested program number at stand-by position by pushing value increase and decrease keys. All values related to the selected program will be shown on the relevant displays. Check that all values are correct.

- Load the rotor with samples by paying attention to the dynamic and static balances.
- Close the lid, see that the lid open warning led (11) turns off.
- Push START key (6) and see that the start led (10) turns on.
- The rotor starts accelerating according to the set value and the acceleration led turns on. The set speed or RCF value is displayed during the set time duration. (RPM or RCF screen will appear in every 3 seconds). For the operations whose speeds are higher than 2.000 rpm, the rotor accelerates according to the set acceleration value up to 2.000 rpm and then continues accelerating with the highest acceleration value, 9 to reach the set speed. As the speed reaches the set value, “-“ sign appears on the acc/brk/temp display (3). (The NF 800 R shows the bowl temperature on the acc/brk/temp display until the centrifugation ends)
- The elapsed time is counted down as the centrifugation starts.
- When the time display shows “00”, the program ends, the breaking led turns on and the set break rate is shown at acc/brk/temp display (3). The rotor brakes with the highest break rate, 9 until it slows down to 2.000 rpm and continues braking according to the set braking value.
- The elapsed time until the rotor stops is counted and is shown at the time display. This is the time that passes until the rotor stops.

- When the speed display shows "0000", the start led (10) turns off, the "end" expression appears and the lid led (9) turns on. The user is warned by an intermittently sounding alarm.
- Push the LID button (8) to open the lid.

You may leave the centrifuge at stand-by position.

**NOTE 1 :** You may display the RPM or RCF value on the speed display during the run by pushing the function button.

**NOTE 2 :** You may see all set values during the run by pushing the value increase key.

**NOTE 3 :** At hold position, the rotor keeps spinning until the STOP key (7) is pushed. After the STOP has been pushed, the rotor brakes according to the set brake rate,

**NOTE 4 :** The centrifuge waits at the stand-by position with the program values of the last program.

**NOTE 5 :** Please open the lid after every operation, otherwise the centrifuge cannot be run again.

## 5.5 SAFETY INTERLOCK SYSTEM

The safety interlock system prevents opening of the lid when the rotor is spinning.

The centrifuge does not operate until the lid is closed and the lid remains locked until the rotor stops spinning. The "Lid Open" indicator turns on to warn the user if the lid is not closed properly.

**NOTE:** If power failure occurs, access to the samples is possible by opening the lid with a special tool. Please see the manual lid opening section (3.6.) for further information.

## 5.6 IMBALANCE DETECTION SYSTEM

- The imbalance detection system of the NF 800 / NF 800 R centrifuges operates electronically when an unacceptable imbalance occurs. In this case the brake is applied immediately.
- The "Err 1" imbalance error is displayed and the imbalance indicator turns on. The lid can not be opened until the rotor stops spinning. The centrifuge can only be started after the lid has been opened and the rotor has been re-loaded correctly.
- To avoid facing imbalance problem, please make sure to insert the tubes correctly.

## 5.7 MANUAL LID OPENING

In case of power cut or any breakdown, the centrifuge can be opened manually to access the samples.

To open the lid manually :

- Switch off the instrument.
- Insert the manual lid opening tool into the hole on the left side of the instrument (see the general view drawing).
- Push the tool keeping it horizontal until the lid opens.

### WARNING !!!

Before opening the lid manually, make sure that the rotor already stopped spinning. Upon opening the lid, lift it by hand and observe the rotor. If the rotor is still spinning, close the lid and wait approximately 10 minutes before repeating the operation. This operation must be carried by someone who is informed of the danger and of the precautions which must be undertaken.

## SECTION 6

### OPERATING PRINCIPLES

#### 6.1 PREPARATION OF THE ROTOR TO RUN

- Before installation, check the rotor for corrosion and cleanliness.
- Chemical corrosion or mechanical corrosion may do severe damage to the rotor and the centrifuge. Particles which are stuck inside the inserts cause the breakage of tubes and lead to major imbalance please check to make sure that no particles are left on the rotor.
- The central hole of the rotor and the motor shaft should also be clean and dry before all centrifugal operations and they must be kept in that way all the time.

#### 6.2 LOADING

The most important condition of an efficient centrifugation is to balance the tubes properly. Loading must be done by meeting the requisites of static and dynamic balance.

**Static Balance:** This balance states that diametrically oppositely replaced weights are almost the same. In application, the liquid level in the tubes should be at the same height to balance the load.

**Dynamic Balance:** This balance states that diametrically oppositely replaced tubes' center of gravity 's are symmetrical with respect to spinning axis of the rotor.

- Although in hospitals the samples having almost the same densities are centrifuged, in industry, samples having different densities may be centrifuged. In this case, the dynamic balance becomes more important factor than the static balance is.
- If the number of tubes to be centrifuged is less than the capacity of the rotor, the tubes must be placed oppositely. If an odd number of tubes is centrifuged, a water filled tube at the same weight should be used for balancing.



**Imbalance of the rotor may cause major damage to the rotor and centrifuge.**



**Never attempt to introduce liquids into the tube inserts.**




**Balance the rotor with glass tubes if you use glass tubes for centrifugation.  
Balance the rotor with plastic tubes if you use plastic tubes for centrifugation.**



**Always use tubes which can withstand to the set speeds.**

Examples of the proper and improper loading are shown below.



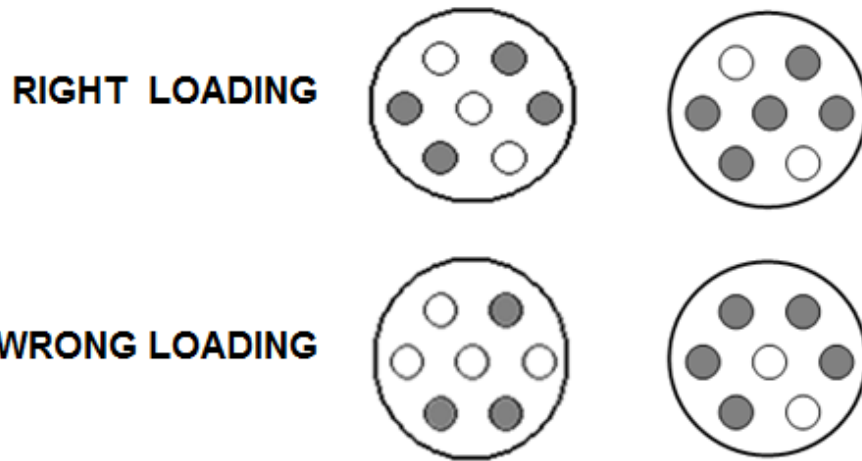


**Figure 3**

**NOTE:** Distribution of the load is more important than the load itself.

### 6.3 ROTOR INSTALLATION

- Put some light oil on the drive head to prevent sticking.
- Be careful that the rotor fits on the drive shaft (See Figure 4).
- Screw the shaft nut with socket wrench to the clockwise direction. Make sure that shaft nut is screwed tightly, but do not expose over-force to the socket wrench while screwing.
- For swing-out rotors, lubricate the pins where the buckets are hanged.
- Place the buckets to the rotor.
- Connect the centrifuge to the power supply.

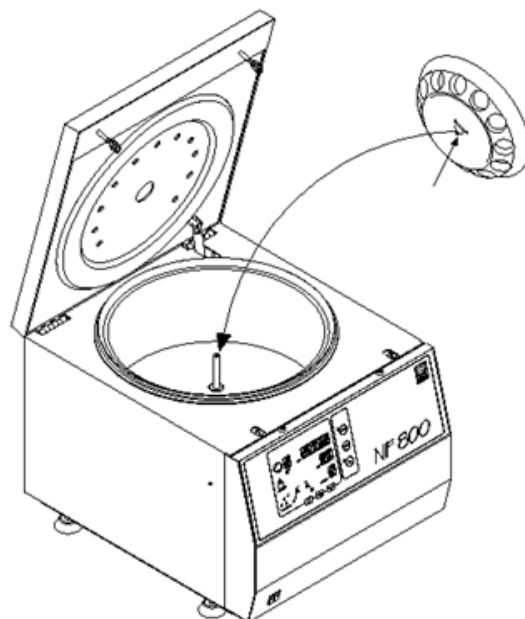

**Figure 4**

## 6.4 DRIVE SYSTEM

- The rotor is driven by a three phase asynchronous motor. The microprocessor control system assures the correct drive speed.
- The force applied to the rotor is directly related to the shape of the rotor, the swing-out rotor receives more load than the angle rotor does. Longer radius and more accessories increase the load of the rotor and decrease the spinning speed.
- The centrifuge does not allow the rotor to spin at a speed which it can not resist mechanically.
- Please make that sure the correct type of rotor is selected during programming.

# SECTION 7

## CLEANING AND PERIODIC MAINTENANCE

### 7.1 PERIODIC MAINTENANCE

- Disconnect the power cable before each periodical maintenance and after the maintenance check the presence of the mains grounding line.
- Rotors should be washed after every use, especially if a spillage has occurred, in warm water containing a few drops of liquid soap. (A mild washing liquid is ideal as a cleaner).
- Rotors and other accessories must be cleaned if any spillage, specially chemicals, occurs
- You may use a nylon brush to clean the buckets and tube inserts of the rotor.
- **Do not use metal brushes.**
- Dry the rotor with a piece of soft absorber cloth. Please make sure that the buckets and inserts are well dried, you may use hair dryer.
- The buckets of the swing-out rotors should be greased frequently with the oil provided with the centrifuge. Please remove the light oil from the pins and put a small amount of fresh oil every time you grease. This will ensure free swinging of the buckets. Most of the imbalance problems are mostly raised by the users who do not clean and oil the pins.
- Please do not leave the rotor on a metal surface, particularly stainless steel as electrochemical reactions set off easily with the aluminum or magnesium in the rotor.
- Make sure that no deposit remains at the bottom of the bucket because the pressure of a flask or tube from above during centrifugation will certainly increase the chance of corrosion.

### 7.2 STERILIZATION

- Apply alcohol , for example %70 ethanol or isopropanol, for 10 minutes against bacteria and viruses.
- The rotors and buckets may be autoclaved at 121°C and under 215 kPa pressure for 20 minutes but please do not forget to remove all accessories.
- Do not use for the sterilization process of formaldehyde.
- Phenol is a corrosive substance and should never be used.
- Glutaraldehyde is a toxic substance and increases the rate of fatty acid in the body.

### 7.3 CORROSION INFORMATION

- Nuve rotors which are made of aluminum are designed to spin at proper RCFs for many years. When used properly, their resistance to corrosion and their life span increases and the imbalance problems decrease.
- All accessories should be checked thoroughly and regularly as almost all laboratories already have the conditions which lead to corrosion easily.

### 7.3.1 Chemical Corrosion

This type of corrosion is caused by chemical reactions. The electrolyte liquid on the surface of the material is the main cause of the chemical reaction. If that electrolyte liquid is allowed to stay at the surface, corrosion occurs. First, discoloration appears and then the metal pitting. Aluminum easily reacts with the ionic solutions.

The other causes of corrosion are as follows,

- Chemical vapors in the laboratory environment which dissolve in the water on the rotor (in refrigerated centrifuges)
- Corrosive liquids which overflow from overfilled and unsealed tubes. (the liquids which spread out during centrifugation)
- Contaminated and non-cleaned buckets, tubes and bottles.

**NOTE:** If the centrifuged samples are corrosive, only rinsing with water is not sufficient. The residuals dissolve in the water and humidity on the rotor and in the buckets.

**NOTE:** Some particles may stick to the tubes, buckets and adapters. These particles crash and do harm to anodized surface during centrifugation and ease the occurrence of corrosion.

### 7.3.2 Stress Corrosion

This type of corrosion is caused by the force of the centrifugation of the corrosive chemical which is already in contact with the alloy. As the aluminum alloy contacts with the corrosive chemical, the stress corrosion starts. This type of corrosion is even more dangerous than the chemical corrosion as the effects of this corrosion are microscopic and very difficult to observe in the course of time.

The corrosive material is pushed against the aluminum alloy by the centrifugation “g” force during the centrifugation. This situation causes the stress corrosion to occur more quickly than the chemical corrosion does. Microscopic cracks occur under the force of the centrifugation.

Every centrifugation causes the aluminum rotor to be attacked by the chemical more and more and eventually micro-cracks decrease the resistance of the rotor against the centrifugation force. Fortunately, no crash occurs just after the first micro-cracks have appeared as the rotors are manufactured according to the high safety limits.

The corrosion of the small amount of corrosive materials does not result in severe cracks but weakens the mechanical resistance of the rotor in the course of time.

The rotor’s places, buckets, buckets’ edges and the base of the rotor should be checked regularly. If needed, they must not be used until a specialist has checked them thoroughly.

## 7.4 CLEANING

- Disconnect the centrifuge before cleaning.
- There is no need of daily cleaning unless a tube breakage occurs or any liquid spills.

## 7.5 ELECTRICITY

Centrifuged at high voltages are present behind the panels. These panels are electrically disconnecting Do not open the centrifuge

## SECTION 8

### DISPOSAL MANAGEMENT CONCEPT

The currently valid local regulations governing disposal must be observed. It is in the responsibility of the user to arrange proper disposal of the individual components.

Applicable local regulations for disposal have to be carefully observed.

The instruments and electronic accessories (without batteries, power packs etc.) must be disposed off according to the regulations for the disposal of electronic components.

Batteries, power packs and similar power source have to be dismantled from electric/electronic parts and disposed off in accordance with applicable local regulations.

## SECTION 9

### TROUBLESHOOTING

**If the centrifuge fails to operate,**

Check that,

- The on/off switch is on,
- The fuses are sound,
- The plug is not defective,
- The centrifuge is well connected to the supply,
- The electricity installation is not defective,
- Power is supplied.

In case of below written failures, related error codes are shown,

#### **Err 1**

- Imbalance. The load is not distributed properly. Balance the load statically and dynamically. Check if the weights of the buckets are the same.

#### **Err 2**

- Motor overheat. Please switch off the centrifuge for at least 20 minutes

#### **Err 3**

- Speed failure. The optical sensor could have failed. Please observe whether the “running” indicator leds are on or not. If it is not on, the optical sensor is defective. If it is on, the motor could have failed. Contact to your service if either occurs.
- This error code occurs when the mains voltage drops below the lower limit.

#### **Err 4**

- Lid failure. This failure occurs when lid is opened during the centrifugation.

**Err 5**

- The temperature sensor endings are broken. This error code appears only on NF 800R.

**Err 6**


- Motor failure. The motor driver unit (inverter) is defective. Switch off the centrifuge and switch on again after the motor has stopped. If you still see the error code, contact to the service.

**PLEASE CONTACT TO AN AUTHORIZED NUVE AGENT TO SEEK TECHNICAL HELP IF AN ERROR OCCURS.**


# SECTION 10

## ELECTRICAL CIRCUIT DIAGRAMS

### 10.1 ELECTRICAL CIRCUIT DIAGRAM OF NF 800

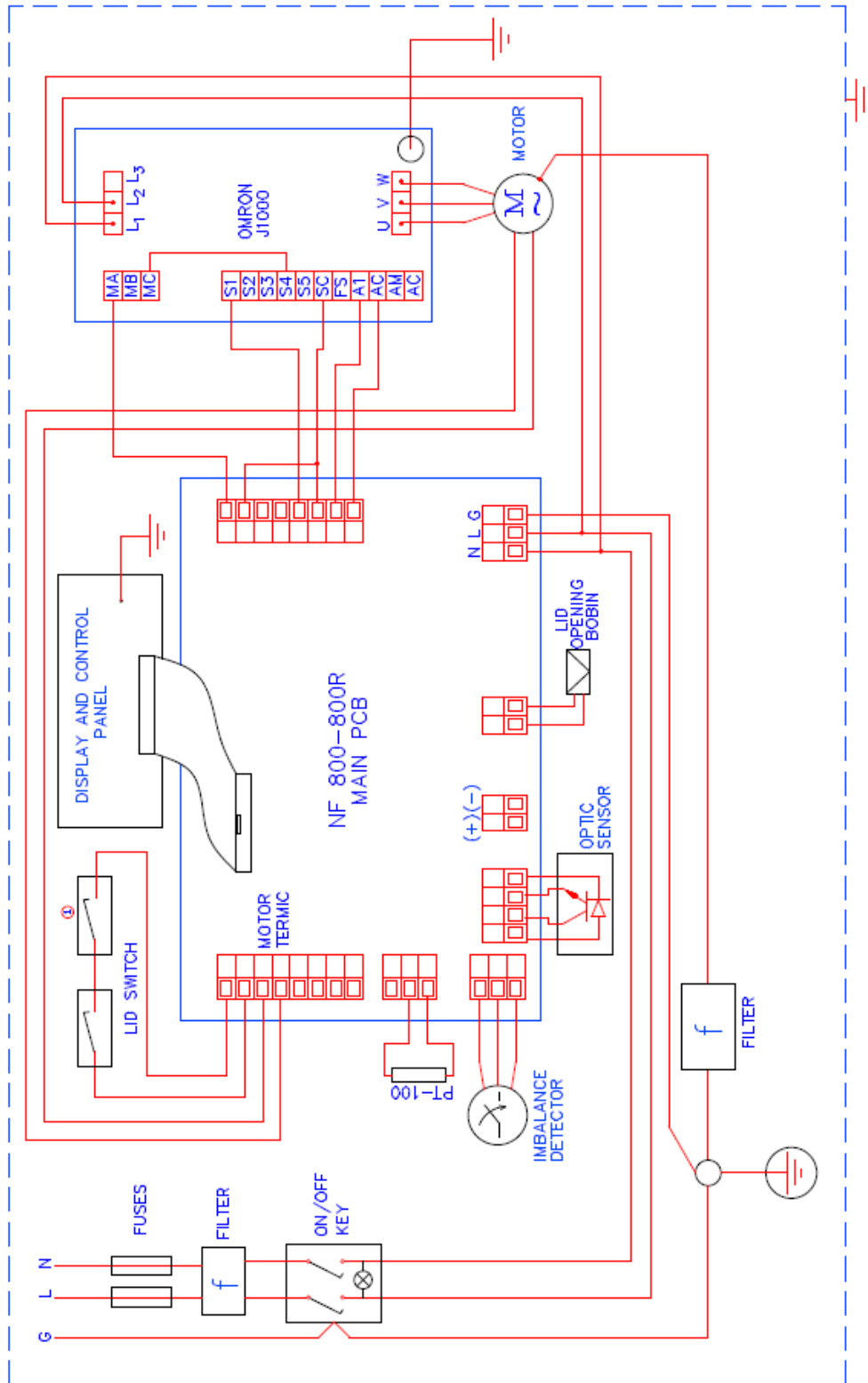

## 10.2 ELECTRICAL CIRCUIT DIAGRAM OF NF 800 R

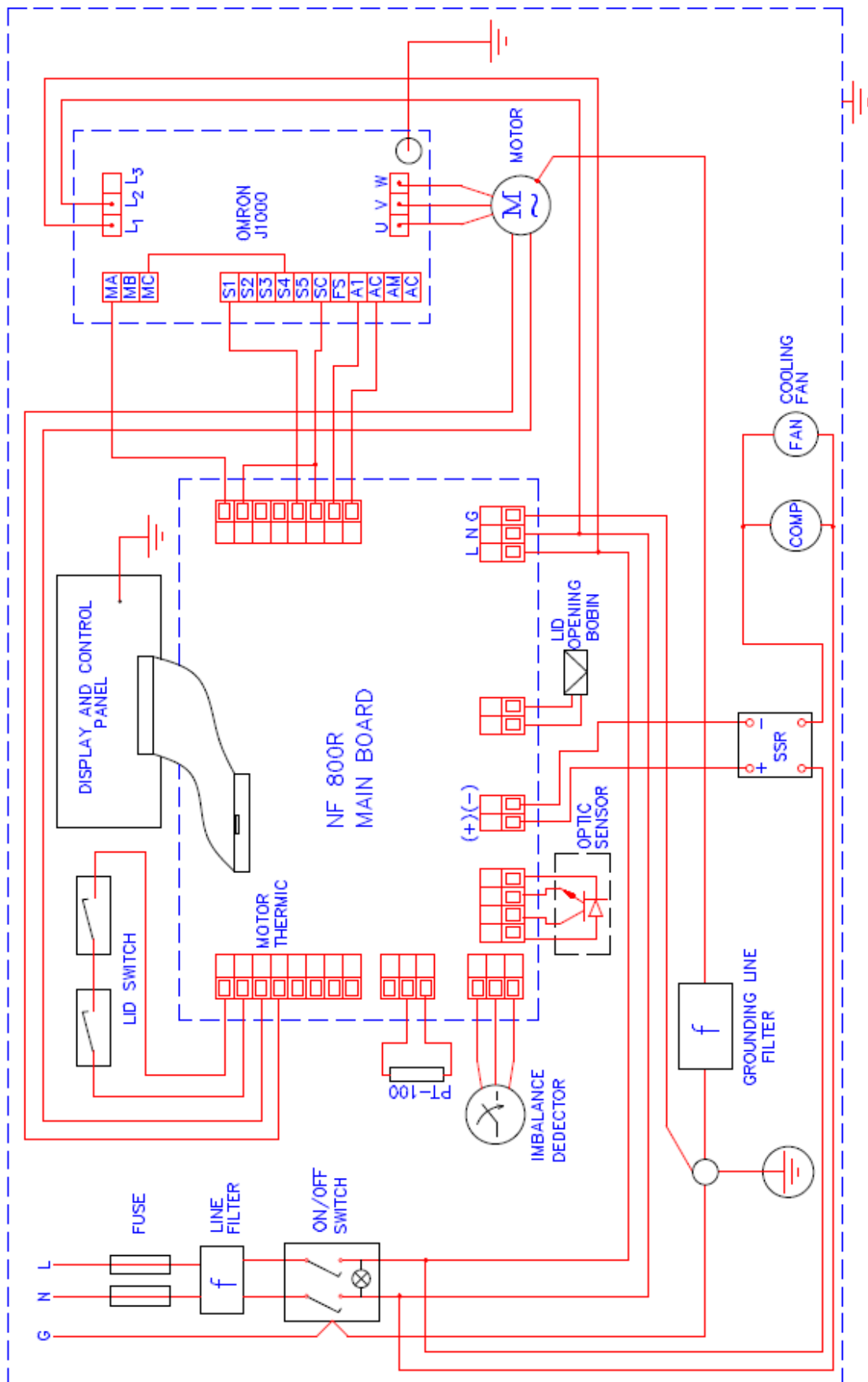

## SECTION 11 WARNING LABEL



GROUNDING PLUG



NF 800R  
FUSES ( 2x10A)



NF 800  
FUSES ( 2x3A)

