



# MICRO CENTRIFUGE TUBES

**S** MALL  
AFE  
IGNIFICANT

Genaxy High-Capacity Micro-centrifuge tubes fit all popular micro-centrifuge rotors. A wide range of coloured options are also available for easy sample identification. The tubes comes in Premium Self standing, Dispenser Packs.

**“1.5ml Micro centrifuge Tubes Are Available With Dual Locking Cap”**

## Unique Features

- » The tubes are manufactured of medical grade polypropylene
- » Manufactured Under ISO Class 8 Clean room Facility
- » Tubes offers exceptional good strength, chemical resistance and withstand temperature from -80°C to +121°C
- » Maximum RCF goes up to 27,000xg with proper rotor support
- » Positive ClickSeal design allows repeated cap closure and eliminates pressure-induced popping of caps during incubation
- » Pierceable cap centre allows easy access to sample by syringe needle.
- » Frosted Writing surface on flat cap and side wall for quick sample identification
- » Superior clarity and external graduations ensure accurate, reproducible lab results
- » Ultra-smooth conical bottoms, which allows for 100% visibility and no sharp “gates” to cause cuts and scratches.
- » DNase/ RNase, Pyrogen, Human DNA, Heavy Metal & Phthalate Free
- » USFDA 21 CFR Compliant

- » Dual locking cap locks in the front of the cap
- » The tube have centrifugation stability upto 30,000 RCF
- » Boil proof, frosted area on exterior & Cap for easy marking




# Genaxy Microcentrifuge Tubes

## Superior Clarity With Greater Strength



Snap Cap  
Long Cap Flanges  
Thin Section For Easy Piercing By Needle  
Flat Cap For Labeling Which Can Be Close With One Hand

Frosted window for writing

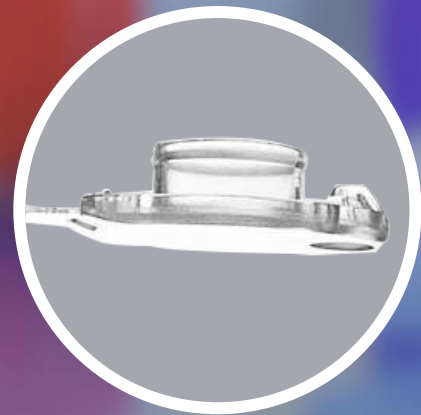

Dual locking cap locks in the front of the cap

Extremely accurate graduation marks

Ultra smooth tube bottoms for 100% visibility - No Sharp Gates to cause cuts and scratches



FOR MORE PRODUCT INFORMATION  
SCAN THE QR CODE



## Micro Centrifuge Tubes

Cat. Number	Description	Qty. / Pack
GEN-MT-060-C	0.6ml Micro Centrifuge Tubes, Clear, Bulk, Non-Sterile	1000 tubes
GEN-MT-060-A	0.6ml Micro Centrifuge Tubes, Assorted, Bulk, Non-Sterile	1000 tubes
GEN-MT-060-Color	0.6ml Micro Centrifuge Tubes, Color*, Bulk, Non-Sterile	1000 tubes
GEN-MT-060-X	0.6ml Micro Centrifuge Tubes, Amber, Bulk, Non-Sterile	1000 tubes
GEN-MT-060-L-C	0.6ml Micro Centrifuge Tubes, Clear, Bulk, Low Retention, Non-Sterile	1000 tubes
*GEN-MT-150-C	1.5ml Micro Centrifuge Tubes, Clear, Bulk, Non-Sterile	500 tubes
*GEN-MT-150-A	1.5ml Micro Centrifuge Tubes, Assorted, Bulk, Non-Sterile	2 X 500 tubes
*GEN-MT-150-Color	1.5ml Micro Centrifuge Tubes, Color*, Bulk, Non-Sterile	2 X 500 tubes
*GEN-MT-150-X	1.5ml Micro Centrifuge Tubes, Amber, Bulk, Non-Sterile	2 X 500 tubes
*GEN-MT-150-L-C	1.5ml Micro Centrifuge Tubes, Clear, Bulk, Low Retention, Non-Sterile	500 tubes
GEN-MT-175-C	1.75ml Micro Centrifuge Tubes, Clear, Bulk, Non-Sterile	500 tubes
GEN-MT-175-A	1.75ml Micro Centrifuge Tubes, Assorted, Bulk, Non-Sterile	2 X 500 tubes
GEN-MT-175-Color	1.75ml Micro Centrifuge Tubes, Color*, Bulk, Non-Sterile	2 X 500 tubes
GEN-MT-175-X	1.75ml Micro Centrifuge Tubes, Amber, Bulk, Non-Sterile	2 X 500 tubes
GEN-MT-175-L-C	1.75ml Micro Centrifuge Tubes, Clear, Bulk, Low Retention, Non-Sterile	500 tubes
GEN-MT-200-C	2.0ml Micro Centrifuge Tubes, Clear, Bulk, Non-Sterile	500 tubes
GEN-MT-200-A	2.0ml Micro Centrifuge Tubes, Assorted, Bulk, Non-Sterile	2 X 500 tubes
GEN-MT-200-Color	2.0ml Micro Centrifuge Tubes, Color*, Bulk, Non-Sterile	2 X 500 tubes
GEN-MT-200-X	2.0ml Micro Centrifuge Tubes, Amber, Bulk, Non-Sterile	2 X 500 tubes
GEN-MT-200-L-C	2.0ml Micro Centrifuge Tubes, Clear, Bulk, Low Retention, Non-Sterile	500 tubes
GEN-MT-060-C-S	0.6ml Micro Centrifuge Tubes, Clear, Bulk, Sterile	100 X 10 tubes
GEN-MT-150-C-S	1.5ml Micro Centrifuge Tubes, Clear, Bulk, Sterile	10 X 50 tubes
GEN-MT-175-C-S	1.75ml Micro Centrifuge Tubes, Clear, Bulk, Sterile	10 X 50 tubes
GEN-MT-200-C-S	2.0ml Micro Centrifuge Tubes, Clear, Bulk, Sterile	10 X 50 tubes
GEN-MT-060-L-C-S	0.6ml Micro Centrifuge Tubes, Clear, Bulk, Low Retention, Sterile	10 X 100 tubes
GEN-MT-150-L-C-S	1.5ml Micro Centrifuge Tubes, Clear, Bulk, Low Retention, Sterile	10 X 50 tubes
GEN-MT-175-L-C-S	1.75ml Micro Centrifuge Tubes, Clear, Bulk, Low Retention, Sterile	10 X 50 tubes
GEN-MT-200-L-C-S	2.0ml Micro Centrifuge Tubes, Clear, Bulk, Low Retention, Sterile	10 X 50 tubes
GEN-MT-5ml-RC	5ml Micro Tube, Natural Bulk Round Bottom, Non-Sterile	1000 tubes
GEN-MT-5ml-VC	5ml Micro Tube, Natural Bulk Conical Bottom, Non-Sterile	200 X 5 tubes

\*GEN-MT-150 - Now Available With Dual Lock System      Color\*\_ Blue, Green, Orange, Red, Voilet, Yellow

### Genaxy Scientific Pvt. Ltd.

155-156, C-Block, Community Centre, Plot 32 to 57, Behind Janak Cinema, Janak Puri, New Delhi-110 058  
 Phones : +91-11-4761 9999 (100 Lines) Fax : +91-11-2553 4163, E-mail : info@genaxy.com