



PETRI DISHES

*TRANSPARENCY
MATTERS*

Genaxy petri dishes are made of Medical Grade virgin polypropylene and is best for microbiological purpose. These plates comes with three ventilation ribs for free air circulation. It also reduces the condensation during incubation. Flared lid skirt and squared corners allows easy one hand operation. Various packaging options are available in the form of single wrap, double wrap, sleeve of 10 and sleeve of 20.

Unique Features


- » Genaxy petri dishes are molded from medical grade virgin polystyrene
 - » Manufactured under ISO Class 8 clean room facility
 - » Ideal for general microbiological use
 - » Clear top & bottom Viewing surfaces
 - » Radiation sterilized : Gamma or ETO
 - » Certificates of Sterilization processing available on request
 - » Available in Various packaging: individual wrap, sleeve of 10 pieces, & sleeve of 20 pieces. Sleeves can be single or double wrap.
 - » Certified DNase/RNase, Human DNA, Phthalate and Pyrogen Free
 - » USFDA 21 CFR Compliant
- Different types of petri plate are available to serve different applications:
- » **90x15mm petri dishes** : comes with Ventilation ribs for free air circulation & reduces condensation during incubation. Consistent flatness allows for even media distribution. Also available with Two & Three compartments for easy identification which is specifically designed for simultaneous studies of organism in a common environment.
 - » **Media saver petri dishes**: comes in smaller diameter thus uses lesser media. Ideal for use in automation. Available in 35mm & 60 mm diameter
 - » **Contact/RHODAC plates**: are useful for bio-burden counts on working surface. Comes with Domed base or flat base with grid molded on the surface of the base. This unique base prevents contamination from fingers touching media. Compatible to various brands of air samplers such as Hi-Media.



Genaxy Petri Dishes



Flared lid skirt for easy one hand operation especially with gloves


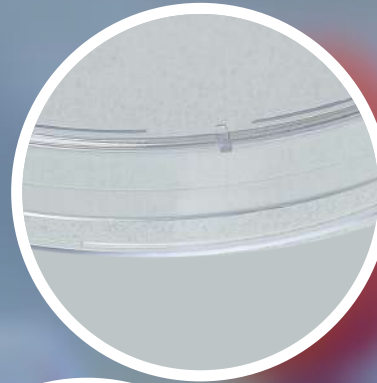


Raised edges on lid for easy stacking



Ventilation ribs allow for free air circulation and reduces condensation during incubation

Consistent flatness allow for even media distribution

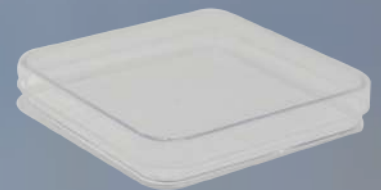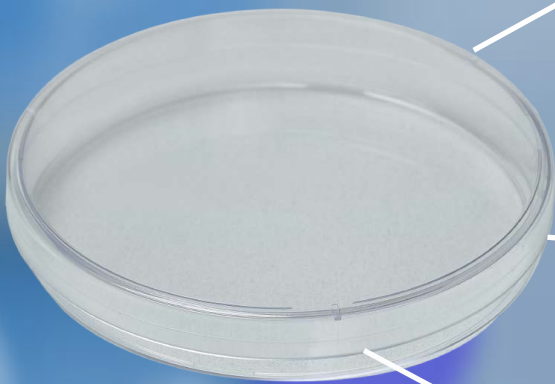

FOR MORE PRODUCT INFORMATION
SCAN THE QR CODE



Petri Dishes

| Cat. Number | Description | Qty. / Pack |
|----------------------|---|-------------|
| GEN-PTD-35 | 35x12mm Petri Dishes, PS, Individual Peel Off, "MEDIA SAVER", Gamma Sterile | 500 |
| GEN-PTD-35-10 | 35x12mm Petri Dishes, PS, Packed in sleeve of 10 Pcs, "MEDIA SAVER", Gamma Sterile | 500 |
| GEN-PTD-60 | 60x15mm Petri Dishes, PS, Individual Peel Off, "MEDIA SAVER", Gamma Sterile | 300 |
| GEN-PTD-90 | 90x15mm Petri Dishes, PS, Individual Peel Off, Fully stackable lid, Gamma Sterile | 460 |
| GEN-PTD-90PP | 90x15mm Petri Dishes, Individual Peel Off, Fully stackable lid, Polypropylene, Autoclavable | 200 |
| GEN-PTD-90-ETO | 90x15mm Petri Dishes, PS, Individual Peel Off, ETO Sterile | 460 |
| GEN-PTD-90-10 | 90x15mm Petri Dishes, PS, Packed in sleeve of 10 Pcs, Gamma Sterile | 540 |
| GEN-PTD-90-10-NV | 90x15mm Petri Dishes, PS, Packed in sleeve of 10 Pcs, Gamma Sterile, Non-Vented | 540 |
| GEN-PTD-90-10-ETO | 90x15mm Petri Dishes, PS, Packed in sleeve of 10 Pcs, ETO Sterile | 540 |
| GEN-PTD-90-10-ETO-NV | 90x15mm Petri Dishes, PS, Packed in sleeve of 10 Pcs, ETO Sterile, Non-Vented | 540 |
| GEN-PTD-90-10-D/W | 90x15mm Petri Dishes, PS, Packed in sleeve of 10 Pcs, Double Wrapped, Gamma Sterile | 540 |
| GEN-PTD-90-10-D/W-NV | 90x15mm Petri Dishes, PS, Packed in sleeve of 10 Pcs, Double Wrapped, Gamma Sterile, Non-Vented | 540 |
| GEN-PTD-90-20 | 90x15mm Petri Dishes, PS, Packed in sleeve of 20 Pcs, Gamma Sterile | 540 |
| GEN-PTD-2C-90 | 90x15mm Petri Dishes, PS, 2 Compartment, Individual Peel Off, Gamma Sterile | 460 |
| GEN-PTD-3C-90 | 90x15mm Petri Dishes, PS, 3 Compartment, Individual Peel Off, Gamma Sterile | 460 |
| GEN-CTP-65 | 65x15mm Contact Dish, Convex Base with Grid | 400 |
| GEN-CTP-65-D | 65x15mm Contact Dish, Convex Dome Base with Grid | 400 |
| GEN-PTD-100-10 | 100x20mm Petri Dishes, PS, Packed in sleeve of 10 Pcs, Gamma Sterile | 300 |
| GEN-PTD-150 | 150x15mm Petri Dishes, PS, Individually Packed, Gamma Sterile | 110 |
| GEN-PTD-130SQ | 130x130x15mm Petri Dishes, PS, Square, Gamma Sterile | 150 |

Genaxy Scientific Pvt. Ltd.