



Stir Bars & Stir Bar Retriever

- Magnetic bar is made from ALNICO, one of the strongest permanent magnets available
- The sintered manufacturing process reduces the chance of decoupling (spin outs), encapsulation in PTFE protects both magnet and what it is stirring
- Chemically inert and easy to clean
- Can withstand temperatures up to 250°C, making them suitable for use with both unheated and heated stirrers
- Low friction octagon shape with a center pivot ring helps create a strong mixing vortex, even with more viscous solutions
- The retrieval bar can recover stir bars even in hot solutions when stirring and mixing is no longer required

10 x 5 mm


20 x 8 mm

25 x 8 mm

28 x 8 mm

30 x 8 mm

35 x 9 mm

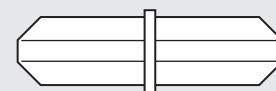

Octagonal Shape

Stir Bars

Pack of 10

Item No.	L x W (in)	L x W (mm)
120547	0.4 x 0.2	10 x 5
120548	0.8 x 0.3	20 x 8
120549	1 x 0.3	25 x 8
120550	1.1 x 0.3	28 x 8
120551	1.2 x 0.3	30 x 8
120552	1.4 x 0.4	35 x 9

Stir Bar Retriever

Pack of 5

Item No.	L (in)	L (cm)
120553	11.8	30

CORPORATE OFFICE
GENAXY SCIENTIFIC PVT. LTD

155-156, C-Block, Community Centre, Plot 32 to 57. Behind Janak Cinema, Janak Puri, New Delhi – 110058
Ph: +91-11-4761 9999(100 Lines) Web: www.genaxy.com E-Mail: technical@genaxy.com