

Multi-purpose, maximum convenience

Innovative Design

- One mixing head does it all with this sleek, modern and compact designed vortex mixer
- Hold slot configurations for multiple sizes of tubes, vessels and assay plates
- Center pad for quick hand-held mixing

Easy to Use

- Control switch with three operation selections (intermittent, off, continuous)
- Variable Speed Rotary Dial provides pre-set speed control selection during shaking mode

Durable Construction

- ABS material resist most commonly used chemicals and is easy to clean
- Stable elastomer feet dampen vibration and prevent movement during operation



Specifications

| | |
|---------------------------|---|
| Speed Range: | To 3,000 rpm, (shaking to vortexing) |
| Orbit: | 3.7 mm |
| Control: | 3 position slider — (top) intermittent, (middle) off, (bottom) continuous |
| Weight: | 8.66 Lb/3.93 Kg |
| Dimensions: | 171w x 190d x 197h mm (6.67w x 7.41d x 7.68h in.) (with mixing head in place) |
| Five-year warranty | |





6 Different Sizes of Tubes
Including standard assay plate



Center pad area for hand held mixing



Assay plate x 1
(to be used with retention cords)



50 mL conical tubes x 2
(horizontally mounted)



15 mL conical tubes x 2
(horizontally mounted)




5 mL tubes x 2
(horizontally mounted in the same retention space as 15 mL)



1.5/2 mL microfuge tubes x 4
(vertically mounted)



0.5 mL microfuge tubes x 6
(vertically mounted)

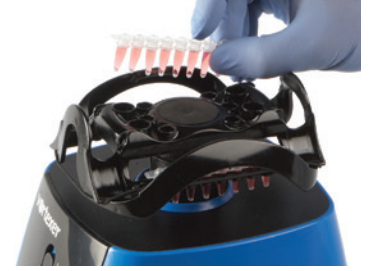

0.2 mL microfuge tubes (strip tubes) x 24
(vertically mounted)

| Item No. | Power Supply | |
|----------|--|--|
| HS120209 | 110V – 120V, 50-60 Hz, 55 Watts, US/NA | |
| HS120212 | 220V – 240V, 50-60 Hz, 44 Watts, CE | |
| HS120214 | 220V – 240V, 50-60 Hz, 44 Watts, UK | |



Corporate Office :

GENAXY SCIENTIFIC PVT. LTD.

155-156, C-Block, Community Centre, Plot 32 to 57, Janakpuri, New Delhi - 110 058

Ph. : +91-11-4761 9999 (100Lines) Web. : www.genaxy.com E-mail : technical@genaxy.com