




MICRO CENTRIFUGE TUBES

S MALL
AFE
IGNIFICANT

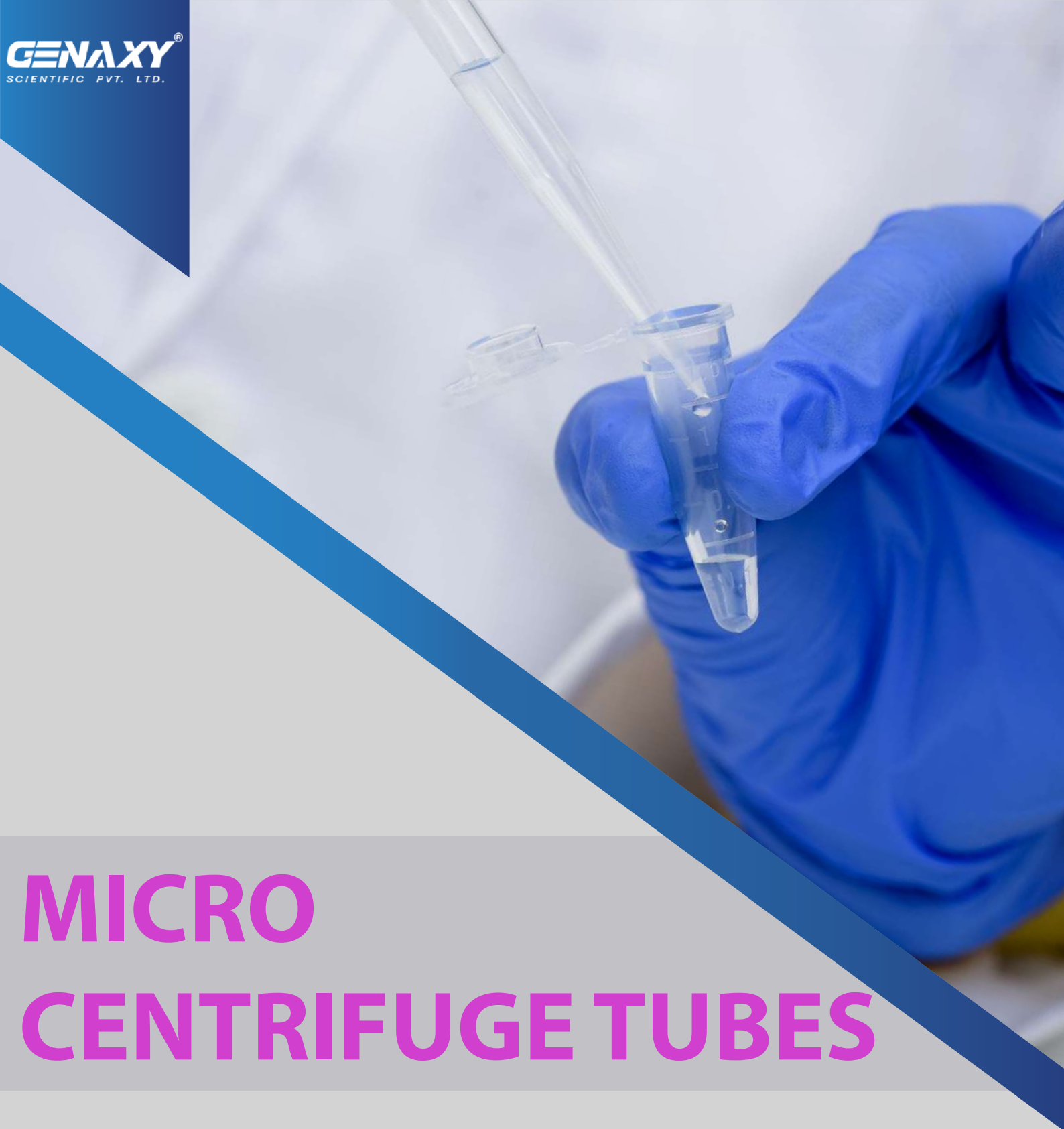

Genaxy High-Capacity Micro-centrifuge tubes fit all popular micro-centrifuge rotors. A wide range of colored options are also available for easy sample identification. The tubes comes in Premium Self standing, Dispenser Packs.

"1.5ml & 2.0ml Micro centrifuge Tubes Are Available With Dual Locking Cap"


Unique Features

- » The tubes are manufactured from medical grade polypropylene which meets USP class IV requirements.
- » Manufactured Under ISO Class 8 Clean room Facility.
- » Tubes offers exceptional good strength, chemical resistance and withstand temperature from -80°C to +121°C. Autoclavable.
- » Maximum RCF goes up to 27,000xg with proper rotor support.
- » Positive ClickSeal design allows repeated cap closure and eliminates pressure-induced popping of caps during incubation.
- » Pierceable cap center allows easy access to sample by syringe needle.
- » Frosted Writing surface on flat cap and side wall for quick sample identification.
- » Superior clarity and external graduations ensure accurate, reproducible lab results.
- » Ultra-smooth conical bottoms, which allows for 100% visibility and no sharp "gates" to cause cuts and scratches.
- » DNase/ RNase, Pyrogen, Human DNA, Heavy Metal & Phthalate Free.
- » USFDA 21 CFR Compliant.
- » Dual locking cap locks in the front of the cap
- » 1.5ml tube have centrifugation stability up to 30,000 RCF and 2.0ml have 27,000 RCF
- » Boil proof, frosted area on exterior & Cap for easy marking



Genaxy Microcentrifuge Tubes

Superior Clarity With Greater Strength

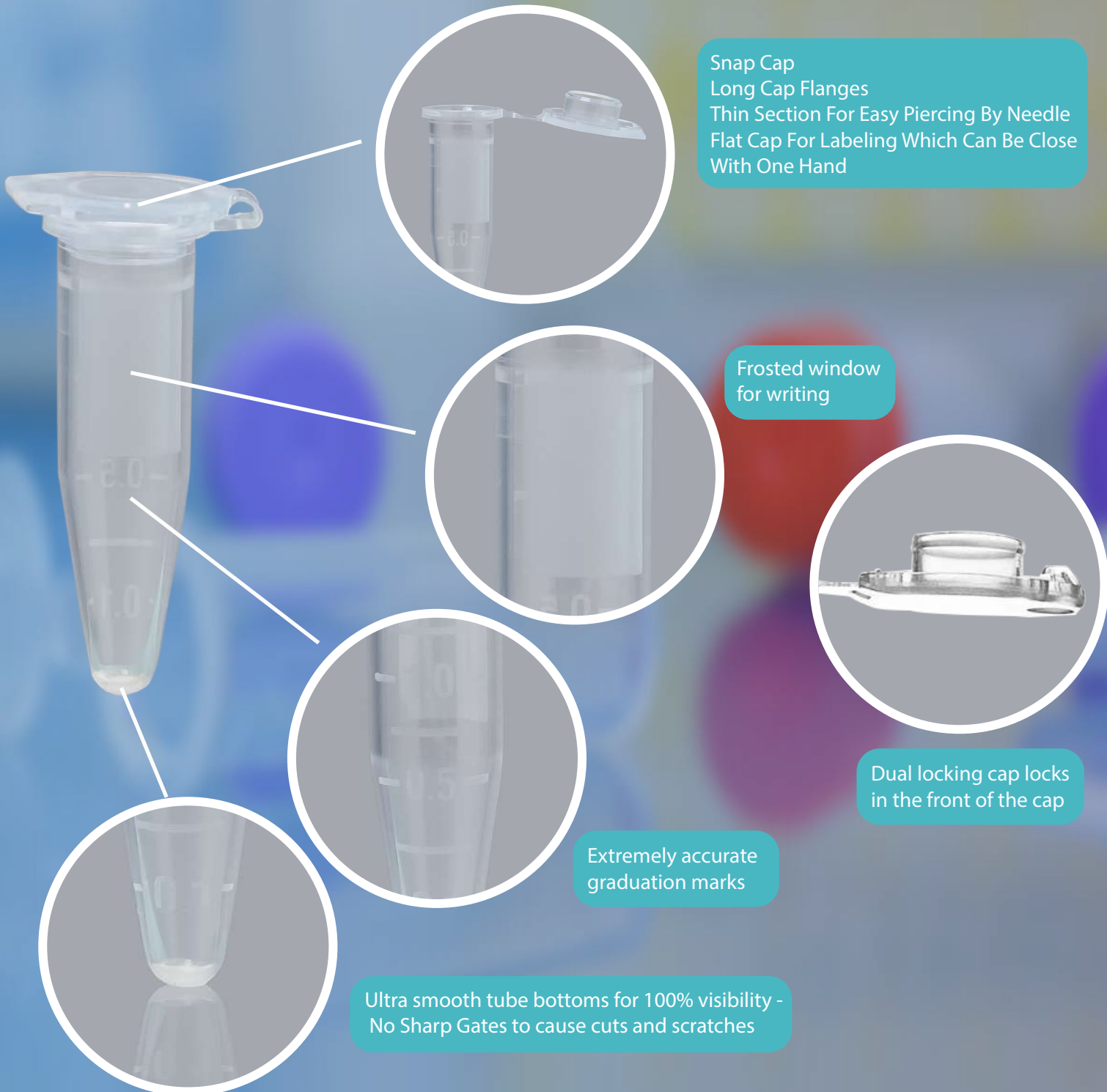

FOR MORE PRODUCT
INFORMATION
SCAN THE QR CODE



Micro Centrifuge Tubes

Cat. Number	Description	Qty. / Pack
GEN-MT-060-C	0.6ml Micro Centrifuge Tubes, Clear, Bulk, Non-Sterile	1000 tubes
GEN-MT-060-C-S	0.6ml Micro Centrifuge Tubes, Clear, Bulk, Sterile	100 X 10 tubes
GEN-MT-060-A	0.6ml Micro Centrifuge Tubes, Assorted, Bulk, Non-Sterile	1000 tubes
GEN-MT-060-Color*	0.6ml Micro Centrifuge Tubes, Color*, Bulk, Non-Sterile	1000 tubes
GEN-MT-060-X	0.6ml Micro Centrifuge Tubes, Amber, Bulk, Non-Sterile	1000 tubes
GEN-MT-060-L-C	0.6ml Micro Centrifuge Tubes, Clear, Bulk, Low Retention, Non-Sterile	1000 tubes
GEN-MT-060-L-C-S	0.6ml Micro Centrifuge Tubes, Clear, Bulk, Low Retention, Sterile	10 X 100 tubes
GEN-MT-150-C	1.5ml Micro Centrifuge Tubes, Clear, Bulk, Non-Sterile	500 tubes
GEN-MT-150-DC	1.5ml Micro Centrifuge Tubes, Clear, Bulk, Non-Sterile, Dual-Lock	500 tubes
GEN-MT-150-C-S	1.5ml Micro Centrifuge Tubes, Clear, Bulk, Sterile	10 X 50 tubes
GEN-MT-150-A	1.5ml Micro Centrifuge Tubes, Assorted, Bulk, Non-Sterile	2 X 500 tubes
GEN-MT-150-Color*	1.5ml Micro Centrifuge Tubes, Color*, Bulk, Non-Sterile	2 X 500 tubes
GEN-MT-150-X	1.5ml Micro Centrifuge Tubes, Amber, Bulk, Non-Sterile	500 tubes
GEN-MT-150-L-C	1.5ml Micro Centrifuge Tubes, Clear, Bulk, Low Retention, Non-Sterile	2 X 500 tubes
GEN-MT-150-L-C-S	1.5ml Micro Centrifuge Tubes, Clear, Bulk, Low Retention, Sterile	10 X 50 tubes
GEN-MT-175-C	1.75ml Micro Centrifuge Tubes, Clear, Bulk, Non-Sterile	500 tubes
GEN-MT-175-C-S	1.75ml Micro Centrifuge Tubes, Clear, Bulk, Sterile	10 X 50 tubes
GEN-MT-175-A	1.75ml Micro Centrifuge Tubes, Assorted, Bulk, Non-Sterile	2 X 500 tubes
GEN-MT-175-Color*	1.75ml Micro Centrifuge Tubes, Color*, Bulk, Non-Sterile	2 X 500 tubes
GEN-MT-175-X	1.75ml Micro Centrifuge Tubes, Amber, Bulk, Non-Sterile	2 X 500 tubes
GEN-MT-175-L-C	1.75ml Micro Centrifuge Tubes, Clear, Bulk, Low Retention, Non-Sterile	500 tubes
GEN-MT-175-L-C-S	1.75ml Micro Centrifuge Tubes, Clear, Bulk, Low Retention, Sterile	10 X 50 tubes
GEN-MT-200-C	2.0ml Micro Centrifuge Tubes, Clear, Bulk, Non-Sterile	500 tubes
GEN-MT-200-DC	2.0ml Micro Centrifuge Tubes, Dual Lock, Clear, Bulk, Non-Sterile	500 tubes
GEN-MT-200-C-S	2.0ml Micro Centrifuge Tubes, Clear, Bulk, Sterile	10 X 50 tubes
GEN-MT-200-A	2.0ml Micro Centrifuge Tubes, Assorted, Bulk, Non-Sterile	2 X 500 tubes
GEN-MT-200-Color*	2.0ml Micro Centrifuge Tubes, Color*, Bulk, Non-Sterile	2 X 500 tubes
GEN-MT-200-X	2.0ml Micro Centrifuge Tubes, Amber, Bulk, Non-Sterile	2 X 500 tubes
GEN-MT-200-L-C	2.0ml Micro Centrifuge Tubes, Clear, Bulk, Low Retention, Non-Sterile	500 tubes
GEN-MT-200-L-C-S	2.0ml Micro Centrifuge Tubes, Clear, Bulk, Low Retention, Sterile	1000 tubes
GEN-MT-5ml-RC	5ml Micro Tube, Natural Bulk Round Bottom, Non-Sterile	1000 tubes
GEN-MT-5ml-VC	5ml Micro Tube, Natural Bulk Conical Bottom, Non-Sterile.	200 tubes
GEN-MT-5ml-VX	5ml Micro Tube, Amber Bulk Conical Bottom, Non-Sterile.	200 tubes

Color* _ Blue, Green, Orange, Red, Voilet, Yellow

Genaxy Scientific Pvt. Ltd.

155-156, C-Block, Community Centre, Plot 32 to 57, Behind Janak Cinema, Janak Puri, New Delhi-110 058

Phones : +91-11-4761 9999 (100 Lines) Fax : +91-11-2553 4163, E-mail : info@genaxy.com